

# Surgical Needles

All needles in this catalog are to scale.



FDA MDR No.: 05148



# **Proper Needle Use**


- DO NOT CLAMP HOLDER DIRECTLY ON THE EYE
- DO NOT GRASP NEEDLE POINT
- DO NOT REUSE BENT OR DAMAGED NEEDLES

Improper surgical technique in the use of a suturing needle may result in needle breakage. Small needle holders should be used on fine gauge needles and heavy holders on heavy gauge needles to avoid undue stress and breakage. As the weakest part of a needle is the eye and the point, Anchor needles are manufactured with a flat area on the shank to facilitate proper clamping and to prevent turning.

Apply holder on flat part about 1/4 of length from the eyed end, do not clamp holder directly on the eye.

When withdrawing needles, clamp on shank about 1/4 length from the point, do not grasp needle point.

Examine needle carefully after each use. If needle becomes bent or damaged during use, withdraw and discard needle.

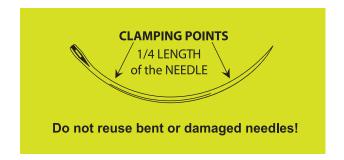



Note: In order to view or print the actual size of needles, set page and printer to A4 and "Do not Scale".

# **SURGICAL NEEDLES SURGICAL NEEDLES TONSIL, DAVIS REGULAR SURGEON'S** 3/8 CIRCLE 1822 1845 1849 2 3/8 Circle 1/2 Circle 4 6 **TROCAR** 1/2 CIRCLE 8 10 1826 12 14 16 CUTTING EDGE 18 20 **REGULAR SURGEON'S** 1/2 CIRCLE 5 1834 TROCAR POINT 6 00 0 3/8 CIRCLE **TROCAR** 1836 2 2 3 TROCAR POINT CUTTING EDGE 4 mc.ph mc.ph


## **SURGICAL NEEDLES**

## **SURGICAL NEEDLES**

# **MURPHY INTESTINAL**

1/2 CIRCLE

1825

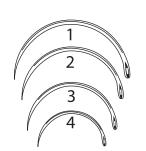



**MAYO INTESTINAL** 3/8 CIRCLE 1837 1 3/4" 1 1/4" 1" 3/4"


# **REGULAR SURGEON'S**

1/2 CURVED


non-sterile only

1821



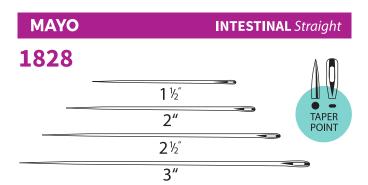


**McGOWAN** 1/2 CIRCLE Heavy Gauge 1859 2

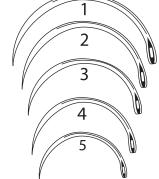




mc.ph

mc.ph

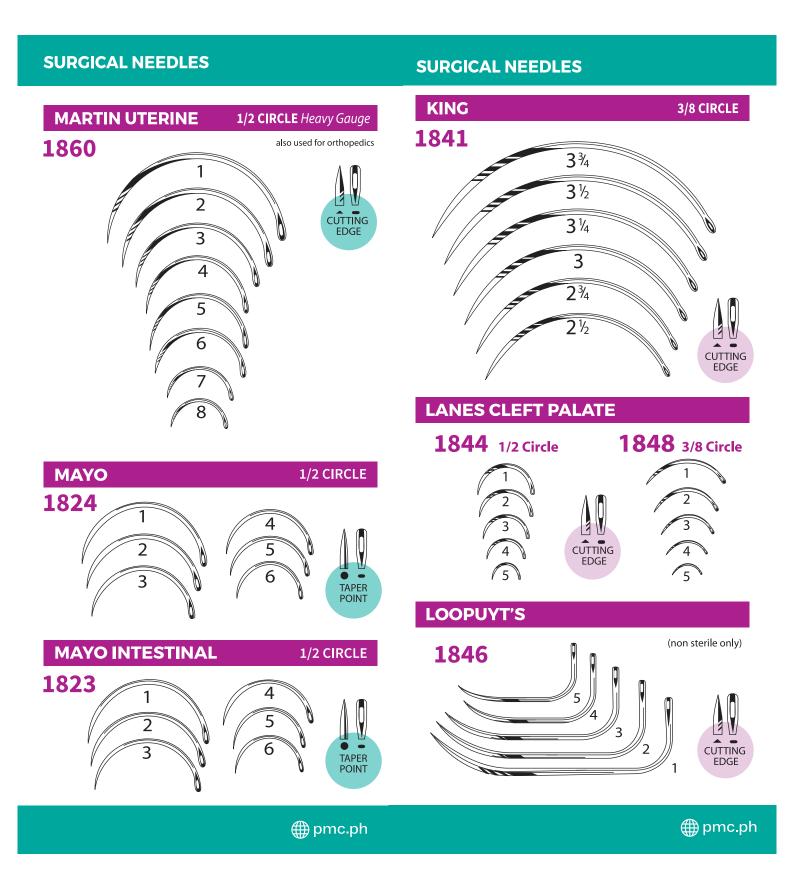

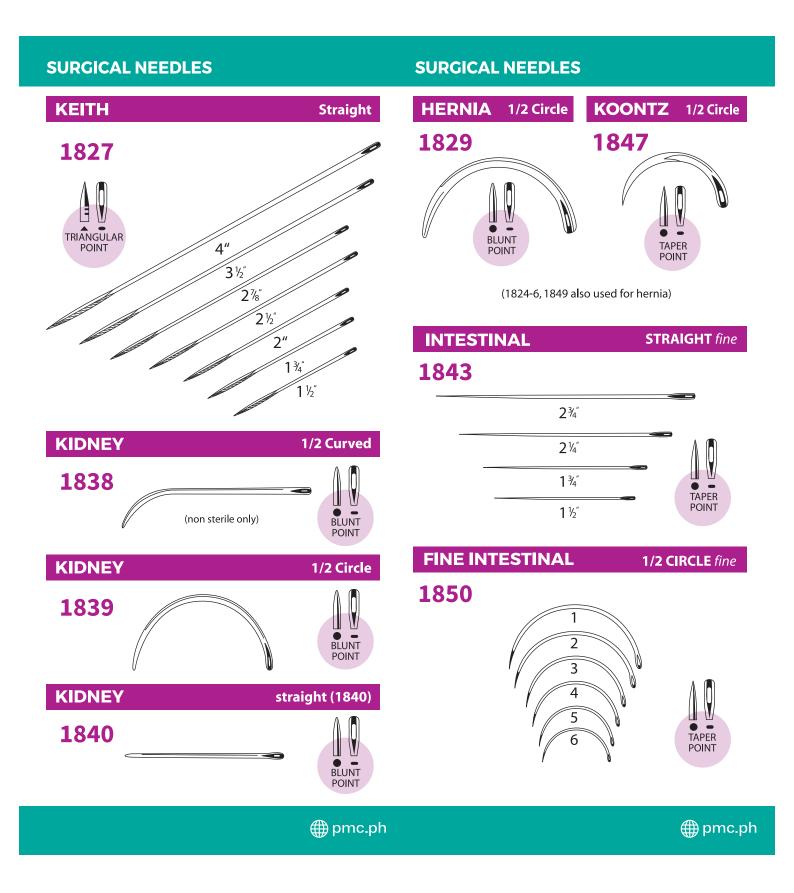

#### **SURGICAL NEEDLES SURGICAL NEEDLES FISH HOOK EYE** 1853-3 **1830** 3/8 Circle **1831** 1/2 Circle CUTTING EDGE **FISTULA** 1/2 CIRCLE heavy gauge CUTTING EDGE 1832 00 1857 Jameson Strabismus **1833** 1/2 Circle (1833-3, 1833-4 also used for dura) (non-sterile only) CUTTING **TAPER POINT EDGE** 3 CUTTING **FERGUSON** 1/2 CIRCLE 1842 **FRENCH SPRING EYE** 1/2 CIRCLE 14 00 1861 2 TAPER 18 4 **POINT /** 20 ` 6 **HAGEDORN STRAIGHT** flat body **22** 8 1852 10 3" 12 21/2" TAPER 2" **POINT** CUTTING EDGE 1 1/2" (1842-20, 1842-22 also used for dura)

⊕ pmc.ph


mc.ph

#### **SURGICAL NEEDLES**

#### **ABSCISSION**

3/8 CIRCLE

## 1851

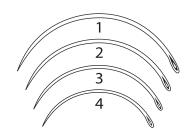




#### **BUNNELS**

**STRAIGHT** 


# 1854

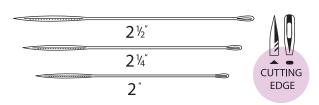

# **CERVIX**

3/8 CIRCLE

# 1835

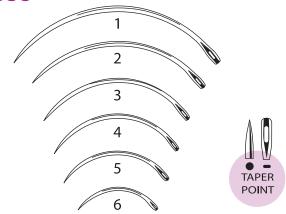

# **Economy and Savings**

TRUE ECONOMY LIES IN QUALITY. ANCHOR STAINLESS STEEL NEEDLES ASSURES BOTH, QUALITY AND ECONOMY.

Discriminating surgeons prefer ANCHOR needles because they satisfy the most exacting requirements of their profession.

Hospitals and institutions discovered the great economy in procuring ANCHOR needles, saving expensive suture material and reducing the quantity of needles purchased in as much as ANCHOR needles can be re-sterilized and reused innumerable times.

ANCHOR needles offers a wide range of impressive features:

- Highest quality stainless steel exposed to sophisticated heat treatment; needles are flexible, but tough and sturdy.
- Each ANCHOR needles is individually ground and inspected to assure a keen cutting edge and point and a large, smooth, burr-free eye. This permits easy threading of suture and prevents suture fraying.
- Each ANCHOR needle is electropolished to assure a point of maximum sharpness allowing greatest ease of penetration and smoothness of passage through tissue.
- ANCHOR needles retain keen cutting edges and points, remain clean and bright inspite of frequent sterilization - will not rust or become discolored.
- ANCHOR needles are well-balanced, with uniformly flat, strong eyes; will not turn in needle-holder.
- The reverse cutting edge type ANCHOR needles have the flat surface of the triangular shape inside of the needle-curve, to minimize trauma.

ANCHOR PRODUCTS COMPANY takes pride in its worldwide reputation as a leading manufacturer of high-quality stainless steel surgical needles for 60 years. Anchor needles have long proven their excellence in thousands of hospitals in the United States and abroad, and are the preferred choice of highly trained surgical-care teams. In craftsmanship, as in surgery, techniques can make the difference; techniques that only human hands can provide. Anchor needles must meet the most exacting and precise product-specifications under strict supervision and verification, complying with rigid quality control procedures. In addition to surgical needles, Anchor offers a wide range of other O.R. related products, such as:

- Eyeless Surgical Needles (atraumatic needles) with sutures attached, in sterilized packaging.
- Custom Needle Kits.
- Reusable stainless steel Sterilizing Rack and Box.





